

PORTABLE WHEELCHAIR SCALES WITH HANDRAIL

DETECTO MODELS 6560 / 6570

- Elegant, curved handrail for patient stability while weighing standing
- May be used for standing or wheelchair weighing (6570 includes seat)
- Column folds down for easy transport
- Integral wheels ease transport
- Ramps rest on floor for easier patient access

*DETECTO
model 6570
with fold-
down seat*



*DETECTO model
6560 (without
fold-down seat)*



SEATED, STANDING, AND WHEELCHAIR WEIGHING ALL IN ONE SCALE



Model
6560

The 6560 and 6570's dual-access ramps go all the way to the ground and allow maximum convenience for loading and unloading patients from either side of the scale. These unique ramps require only half the effort for a healthcare professional to roll a wheelchair-bound patient onto the scale than competitor's models.

- Clinical-grade accuracy to 0.2 lb / 0.1 kg
- Quick and easy mobility for transport or storage
- Up to 99 tares (optional) to store wheelchair weights for efficient patient weighing
- Hold/Release key to hold weight after leaving the scale for better patient care
- Interface with a PC or printer via serial, USB or optional Wi-Fi or wired Ethernet for EMR
- Patient ID barcode scanning sends 14-digit patient ID to indicator for EMR/EHR records



Curved Handrail for Patient Stability

The 6560 and 6570's elegant curved handrail provides patient support while weighing standing on the platform.




Single-Pass Weighing Saves Time

Enjoy the time-saving benefits of these scales by using the convenient pushbutton or keypad tare to remove the wheelchair weight.



Fold-up Chair Seat

The 6570's chair seat easily folds up when not in use, allowing wheelchair weighing (not available on 6560).

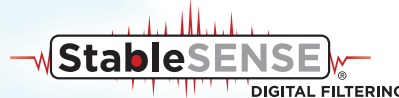

Model
6570

The scale's columns fold down and latch securely for easy transport and vertical fold-up storage requiring minimal floor space.



FOR USE IN **PORTABLE MULTIPLE CLINICAL DEPARTMENTS**

Digital Multifunctional Medical Scales with Time-Saving and Patient Comfort Benefits

DETECTO's 6560 and 6570 portable wheelchair scales offer time-saving and patient comfort benefits. The scales features 1,000 lb / 450 kg capacity, optional Wi-Fi for EMR/EHR, vertical fold-up storage, two-way built-in ramps for accessibility from both sides, large 32 in x 36 in / 81 cm x 91 cm platform, clinical-grade accuracy to 0.2 lb / 0.1 kg, battery or AC power, and up to 99 tares (optional) to store wheelchair weights for efficient patient weighing. Enjoy the time-saving benefits of these scales by using the convenient pushbutton or keypad tare to remove the wheelchair weight. OIML EC class III models are available (see back page).



The scale's columns fold down and latch securely for easy transport and vertical fold-up storage requiring minimal floor space.

EFFICIENT MEASUREMENT



The 6560 and 6570's MedVue® weight indicator is the clinical industry's most advanced, featuring EMR/EHR-ready serial, USB, and optional Wi-Fi connectivity, BMI calculation, yellow "quick keys" for basic weighing, and a multiline LCD display. The scale is powered by 6 C batteries (not included) or optional AC adapter.

Easily-Collapsible Column

By sliding the release lever, you can easily fold down the columns for easy transport or vertical storage.



Low-Profile, Thin Platform

The platform is only 2.2 inches/56 mm high for easy patient accessibility and comfort while weighing either standing or in a wheelchair.



VERSATILE EMR-READY PORTABLE SCALES



6560



6570

	6560	6560-AC	6560-PR-AC	6560KGEU	6570	6570-AC	6570-PR-AC	6570KGEU
Capacity:	1,000 lb x 0.2 lb / 450 kg x 0.1 kg			300 kg x 0.1 kg or 350 kg x 0.2 kg	1,000 lb x 0.2 lb / 450 kg x 0.1 kg			300 kg x 0.1 kg or 350 kg x 0.2 kg
Description:	AC adapter and printer not included	AC adapter included	AC adapter and P600 tape printer with cable included	AC adapter and printer not included	AC adapter and printer not included	AC adapter included	AC adapter and P600 tape printer with cable included	AC adapter and printer not included
OIML EC Class III Approved:				●				●
Platform Size:	32 in x 36 in x 2.2 in / 81 cm x 91 cm x 5.6 cm							
Keypad:	19 key, color-coded, sealed membrane							
Connectivity:	RS232 serial, USB and optional wired and Wi-Fi Ethernet							
LCD Display Dimensions:	3.5 in W x 3 in H / 89 mm W x 76 mm H backlit blue LCD							
Power:	6 C cell Alkaline, Ni-Cad or NiMH batteries (not included) OR 12 VDC 1.25A wall plug-in AC power adapter, Detecto part number MV1PWR (included with -AC models listed above)							
Resolution:	5,000 divisions, max.							
Operating Temperature:	14° to 104°F (-10° to 40°C)							
Ramp Size:	36 in (91.4 cm) W x 7.9 in (20.2 cm) L each							
Display:	5 digit, 7 segment, 7/8 inch (22.23 mm) high LCD weight digits							
Country of Origin:	USA							
UPC Code:	809161302202	809161317008	809161317107	809161304503	809161302301	809161317206	809161317305	809161304602
Net Weight:	116 lb / 53 kg	117 lb / 53 kg	120 lb / 54 kg	116 lb / 53 kg	140 lb / 63 kg	141 lb / 64 kg	144 lb / 65 kg	140 lb / 63 kg
Shipping Weight:	126 lb / 57 kg	127 lb / 58 kg	130 lb / 59 kg	126 lb / 57 kg	150 lb / 68 kg	151 lb / 69 kg	154 lb / 70 kg	150 lb / 68 kg

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



102 E. Daugberty, Webb City, MO 64870 USA
 Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153
 www.Detecto.com

Mobile: m.detecto.com

SOLD BY: