

6745 DIGITAL BABY SCALE



DETECTO model
6745 baby scale

DETECTO®

- 1 gram divisions
- 30 lb / 15 kg capacity
- Roomy 22 in x 9 in / 57 cm x 23 cm weighing tray (flat portion)
- Full-color, touchscreen 3.2 in (diagonal) display with 320 x 240 resolution
- Fun, colorful animal graphics attractive for children
- (6) C batteries (not included) or optional AC adapter
- Integral measuring tape (22 in / 56 cm long)
- Wide supportive scale base for infant safety
- 8-digit numeric patient ID for EMR/EHR
- Memory key to store and recall a weight
- Built-in carrying handles



NTEP Approved




Bluetooth™

Wi-Fi and Bluetooth 5.0 Connectivity



Units locking into lb or kg

6745 | DIGITAL BABY SCALE

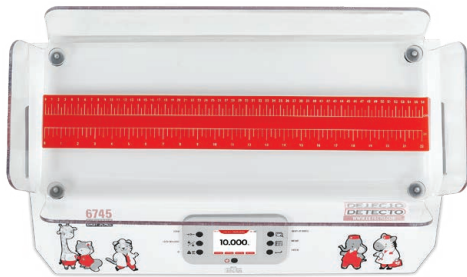

MEASURING TAPE

0-22 in / 0-56 cm length measuring tape imprinted underneath the clear thermoplastic platform, so will never wear off during use.

PATIENT ID				
12345678	1	2	3	
Key in the numeric patient ID to be transmitted when the Print key is pressed.	4	5	6	
	7	8	9	
SAVE	CANCEL	0	.	DEL

PATIENT ID

The 6745's touchscreen graphics display allows you to enter up to an 8-digit numeric identification number for the patient, which can be sent to an EMR/EHR system.



CONNECTIVITY

Ports include 9-pin RS232 for sending data to a serial printer and USB-B for connecting with a PC for EMR/EHR, plus Wi-Fi and Bluetooth wireless connectivity is standard.

Capacity	6745: 0 - 10 lb x 0.002 lb / 10 - 30 lb x 0.005 lb 0 - 10 lb x 0.02 oz / 10 - 30 lb x 0.05 oz 0 - 5 kg x 0.001 kg / 5 - 15 kg x 0.002 kg 0 - 5,000 g x 1 g / 5,000 g - 15,000 g x 2 g	6745KG: 0 - 5 kg x 0.001 kg / 5 - 15 kg x 0.002 kg 0 - 5,000 g x 1 g / 5,000 g - 15,000 g x 2 g
Platform Size (Flat Portion)	22.5 in W x 9.1 in D / 57.1 cm W x 23.1 cm D	
Overall Dimensions	26 in W x 15.1 in D x 7.5 in H / 66 cm W x 38.4 cm D x 19.0 cm H	
Display Type	3.2 in / 81.3 mm (diagonal) full-color TFT display 320 x 240 resolution	
Touch Screen	Resistive touch panel covering the display	
Keys	Mechanical switch type, Power, Zero, Lock/Release, ID, Cycle Display Mode, Menu, Print	
Power Requirements	Six C-size Alkaline, Ni-Cad, or NiMH batteries (not included) Optional 100 to 240V AC 50/60Hz 12V DC 1.2A Medical wall plug-in UL, CSA, TUV, CE Approved AC power adapter (Cardinal part number 6800-1006) and Power Cord, 18GA, detachable IEC 320-C13, UL, CSA Approved, US plug (Cardinal part number 6980-1052)	
Operating Environment	Operated Temperature Range: 50 to 104°F (+10 to +40°C) Humidity: 0 to 90% non-condensing	
Communication Interfaces	Wi-Fi, Bluetooth 5.0, RS232, and USB-B	
Shipping Weight	35 lb / 15.8 kg	



StableSENSE adjustable filtering to compensate for baby movements



SOLD BY:

DETECTO®