

DETECTO®


750

CLINICAL WEIGHT INDICATOR


- Extra-Large 0.75-in/19-mm high LCD Digits
- Units Locking into Pounds or Kilograms Helps Prevent Dosing Errors
- Body Mass Index (BMI) Calculation
- Wi-Fi/Bluetooth Models Available for Wireless EMR/EHR Connectivity
- (1) RS232 Serial Port
- Models are Available with AC Adapter Included
- Lock/Release Key for Viewing Weight After Patient Steps Off the Scale
- Durable Steel Construction for Years of Reliable Service



Wi-Fi/Bluetooth models are available for wireless EMR/EHR connectivity.

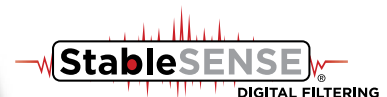

SCALE CONNECTIVITY FOR PATIENT RECORDS

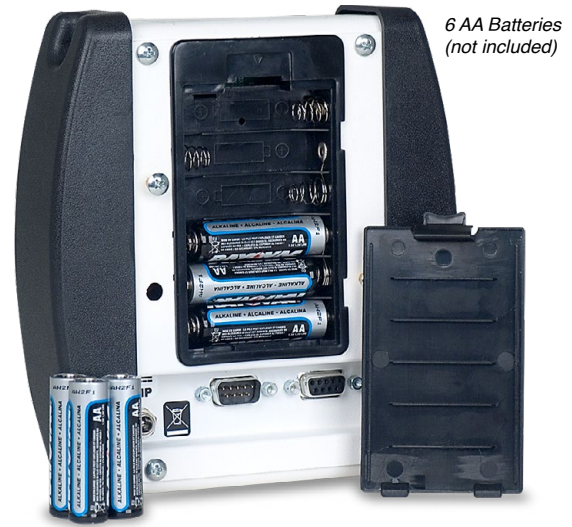

Model	750
Power Requirements	6 AA Alkaline, Ni-Cad or NiMH batteries (not included) OR an optional 100 to 240 VAC 50/60Hz 12 VDC 1A wall plug-in UL/CSA listed AC power adapter (DETECTO part number 6800-1045). Models are available with the AC adapter included.
Operating Temperature	14 to 104°F (-10 to 40°C)
Display	0.7 in / 19 mm high LCD, 6 digit, 7 segment
Sensitivity	0 to 2.4 mV/V
Signal Input Range	0 to 12 mV max
Transducer Excitation	5.0 VDC
Resolution	2,500 divisions
Capacities	Up to 10,000 divisions non-commercial
Division Value	1, 2, or 5
Sample Rate	1 to 16 samples per second
Auto Zero Range	0.5 or 1 through 9 divisions
Weighing Units	Pounds or Kilograms with Units Locking
Keyboard	Membrane type with 5 keys and 4 directional arrows
Enclosure Size	6.125 in W x 6.75 in H x 1.75 in D / 16 cm W x 17 cm H x 4 cm D
Construction	Painted Steel
Connectivity	RS232 serial port standard. Wi-Fi/Bluetooth models are available for wireless EMR/EHR connectivity (model 750C)
Standard Features	Units Locking, Auto Shut-Off, Selectable Sleep Mode, Selectable Filtering, Serial I/O (Input/Output)
UPC Code	809161149609
Country of Origin	USA



DETECTO's model 750 clinical indicator features an easy-to-use keypad and extra-large display making it perfect for patient weighing and viewing by the attending nurse, physician, or patient. A quick-and-easy BMI calculator allows you to view Body Mass Index. The durable steel construction and easy-to-replace AA batteries make the 750 perfect for all clinics, hospitals, and physician's office scales. Wi-Fi/Bluetooth models are available for wireless connectivity with EMR/EHR. Units locking into pounds or kilograms helps prevent dosing errors. Models are available with AC adapter included.



DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com