

DETECTO®

8432-CH SERIES

DIGITAL BABY SCALES WITH INCLINED SHELL SEAT



Infants readily accept the inclined seating of DETECTO's 8432-CH series digital pediatric scales. The built-in, shell-shaped seat is made of plastic material and fun for children to sit in while weighing.

With pound and kilogram readings and heavy-duty steel base and understructure for support and stability, these pediatric scales have a lot to offer.



- Digital speed and accuracy
- 44 lb / 20 kg capacity
- Clinical accuracy to .05 lb / .02 kg
- Simple-to-use weight indicator
- Battery or AC power
- Solid steel base for steady assurance while weighing
- Pushbutton tare removes blanket weights
- Units locking into pounds or kilograms for safe dosages

8432-CH SERIES

DIGITAL BABY SCALES WITH INCLINED SHELL SEAT



Model	8432-CH	8432-CH-AC	8432KG-CH
Capacity	44 lb x 0.05 lb / 20 kg x 0.02 kg	44 lb x 0.05 lb / 20 kg x 0.02 kg	20 kg x 0.02 kg
AC Adapter		•	
Scale Height	26 in / 66 cm		
Display	Six-digit, seven segment, 0.7 inch (17.8 mm) high LCD		
Power	Six "AA" Alkaline, Ni-Cad or NiMH batteries (not included) OR an optional 100 to 240V AC 50/60Hz 12V DC 1A wall plug-in UL/CSA listed AC power adapter (Cardinal part number 6800-1045). Models are available with the AC adapter included.		
Keyboard	Membrane type with 5 keys and 4 directional arrows		
Connectivity	RS232 serial port		
Shipping Weight	33 lb / 15 kg		

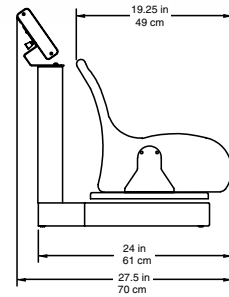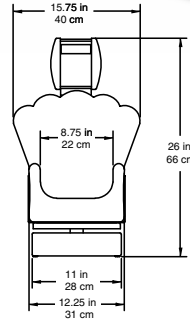

The 8432-CH has a lock-in weight feature which compensates for typical baby movement. The crisp LCD features 0.75 inch high digits for easy readouts.



DETECTO's plastic shell-shaped inclined seat on the 8432-CH series baby scales makes weighing fun for pediatric patients.



The weight indicator is conveniently situated over the top of the shell seat, so it's easy to read the weight while the baby is seated.



DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com