


OPERATING ROOM PRODUCTS

When accuracy is critical, and seconds count, DETECTO scales provide the quick, precise information you need



PR Series

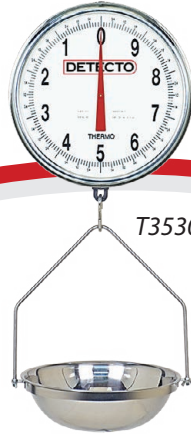

T3530



420-2000



RESCUE
medical carts



GH3SS



WPS12UT

DETECTO®

ACCURATE RELIABLE SCALES AND EQUIPMENT FOR OPERATING ROOMS THAT CAN HANDLE ANYTHING

Dependability is required for equipment in an Operating Room. You know you'll have medical equipment you can count on from DETECTO. The diverse line of hospital scales will meet the demands of the OR, where efficiency and accuracy are critical. Our mobile medical carts, glove boxes, and waste receptacles are built to perfectly serve your needs and perform reliably for a lifetime.

WPS12UT

**Submersible
Digital Scale with
Utility Bowl**



The WPS12UT digital scale with utility bowl is IP67-rated for water protection and can be completely submersed underwater for easy cleaning. This stainless steel scale is ideal for weighing soiled or wet diapers, chux pads, lap sponges, and other materials. The 4.6 in / 11.7 cm deep bowl aids sample retention and allows easy clean-up. The scale's digital speed and accuracy means you'll be able to view weights instantly with up to a fine 0.1 ounce or 1 gram graduations.

420-2000

**Digital Precision
Balance Scale up to
2000 g Capacity,
with 20 mg
Precision**



DETECTO's digital precision balance scales are precise and durable, measuring in a wide variety of units for many applications where accuracy is demanded, such as weighing blood bags, blood and blood components, patient medical waste samples, medication, and other critical materials. The scale has a built-in rechargeable battery and a removeable air shield for precise measurement. The stainless steel platform size is 5.7 in W x 4.5 in D.

T3530

**Hanging Dial Scales with
Stainless Steel Bowl and
32 lb Capacity**



DETECTO's hanging dial scales feature a deep 13.5-in/34-cm diameter removable stainless steel bowl making them ideal for medical weighing use. Protective glass plates cover the 14.5-in/37-cm diameter dials with extra-large readouts visible from both sides. Available in lb or kg models: 32 lb or 15 kg capacities. Full capacity is reached in 3 revolutions of the dial. Thermolastic springs guarantee long-term accuracy while an air dash pot absorbs shock when loads are dropped on the scale protecting the mechanism and increasing the scale life.

RC33669BLU-L

**RESCUE Steel
Medical Carts**



DETECTO's Rescue medical carts are sturdy and reliable... built to make sure you have the supplies you need where and when you need them. They come in several colors, with blue selected for OR use. Carts have your choice of five or six drawers that can be secured with a key or EMG breakaway lock. The ABS countertop is easy to clean, and 5-inch wheels feature two parking locks and one steering lock. Optional accessories include IV poles, cardiac boards, catheter holders, drawer dividers, glove dispensers, sharps container holders and waste bins.

P-100R

**Step-On Cans, the Rugged
Waste Receptacle**



Named for their rugged toughness, Step-On Cans are perfect for hospitals, doctors' offices, patients' rooms, and laboratory settings. Equipped with a heavy-duty foot pedal and linkage and a bag securing mechanism, they come in a variety of five different sizes and capacities (4, 6, 8, 12 and 25 gallon) and two different materials: stainless steel or a baked epoxy finish (available in red or white). The cans are flame retardant, have a self-closing lid for sanitation, and feature unitized construction.

GH3SS

**Stainless Wall-Mount
Glove Box Holders**



Conserve counter space with DETECTO's stainless steel universal glove box holders. The open design allows you to see what size or type of glove is being dispensed. The holders accommodate most glove boxes and are universal enough to be used for tissue holders where counter or shelf space is limited. Holds up to three boxes of gloves. Lifetime warranty.

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:

DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com