

### STRE+CHER

SCALE



### MODELS 8500/8550

- Spacious platform accommodates all common stretcher types
- Two integral wheels and guide handles for mobility
- High-capacity digital scale for stretchers and wheelchairs
- Available with or without column for weight indicator



## PRECISION WEIGHING





- Spacious 60 in W x 32.5 in D / 1,524 mm W x 826 mm D platform size
- Hold/Release key to hold weight after leaving the scale for better patient care
- Sturdy steel construction for optimal durability
- Mobile for use in multiple clinical departments
- Two integrated ramps included for easy access






Detecto

Model 8550

**DETECTO's Stretcher** Scale® is ideal for emergency arrivals into a hospital when time is of the essence for gathering patient weights. The **DETECTO** models 8500 and 8550 can be used in the ER, dialysis, long-term care, and nursing homes.

The single-pass weighing feature saves time by allowing you to enter the stretcher or wheelchair tare weight using the numeric keypad on the indicator. You may also use the convenient pushbutton tare feature to remove the stretcher or wheelchair weight if you don't know the weight already.

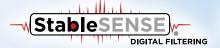







### PORTABLE WITH INTEGRAL WHEELS AND GUIDE HANDLES

DETECTO's Stretcher Scale® is engineered and manufactured at a light enough weight that it is easily portable to lift up and roll by nursing

professionals. The model 8500 is 115 lb/ 52 kg net weight and the model 8550 is 130 lb/59 kg. The slim profile scale design fits easily through doorways of nearly any size.



# STRE+CHER SCALE

### **LIFT AND ROLL**

The scale's relatively low net weight, wheels, and guide handles all combine to make it easy to maneuver.



### **WEIGH WHEELCHAIRS, TOO**

The Stretcher Scale ® may be used to weigh wheelchairs as well as stretchers and gurneys.



### **LOW PROFILE PLATFORM**

DETECTO's Stretcher Scale® has a 2-inch-high, low-profile platform with textured mat surface.



### **BRACKET INCLUDED (Model 8500)**

The model 8500 comes with a wall or desk mount swivel bracket for the MedVue® weight indicator.



### **EASILY-COLLAPSIBLE COLUMN**

It's easy to quickly slide the release lever and fold down the 8550's column for transport.



### **TWO-WAY QUICK ACCESS RAMPS**

The two-way-access ramps rest on the floor making it easy to roll patients onto the scale platform.



### **OPTIONAL GUIDE RAILS**

An optional yellow guide rail kit is available for the 8500/8550 (8500-GUIDERAILS).



SCALE®

DETECTO MODELS 8500/8550








DETECTO model 8550 shown with optional yellow guide rails

### **Options**

8500-GUIDERAILS	Guide rail kit painted yellow for models 8500/8550 (includes wheel side and handle side)
MV1PWR	Medical-grade AC adapter power supply with US plug
MVWIFI	Wireless option card
MVIP	Wired Ethernet interface option

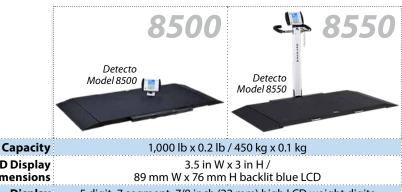

Capacity	1,000 lb x 0.2 lb / 450 kg x 0.1 kg		
LCD Display	3.5 in W x 3 in H /		
Dimensions	89 mm W x 76 mm H backlit blue LCD		
Display	5 digit, 7 segment, 7/8 inch (2	2 mm) high LCD weight digits	
Resolution	5,000 divis	ions, max.	
Keypad	19 key, color-coded, sealed membrane		
Functions	Tare, Net/Gross, Hold/Release, Zero, metric conversion,		
	ID storage, Height entry for Body Mass Index calculation, and printer output		
Power	Optional model MV1PWR AC adapter or 6 C batteries (not included)		
Battery Life	7,560 weighments (with backlight fully ON)		
buttery Life	10,620 weighments (with backlight fully OFF)		
Connectivity	RS232 serial, USB and optional wired and Wi-Fi Ethernet		
Cable Length (Model 8500 only)	6 ft / 1.8 m coiled cable, 14 ft / 4.2 m fully extended (from scale to indicator)		
Mounting Bracket	350-degree indicator swivel mounting bracket included		
Operating	14 to 104 °F /		
Temperature	-10 to +40 °C		
Construction	Mild steel		
Included	Two ramps (assembly required)		
Overall	76 in W x 38 in D x 2 in H /	76 in W x 38 in D x 44.8 in H /	
Dimensions	1,930 mm W x 965 mm D x 51 mm H	1,930 mm W x 965 mm D x 1,138 mm H	
Remote Indicator	9.8 in W x 8.3 in D x 7.8 in H /		
	249 mm W x 211 in D x 198 mm H	Included in	
	(including swivel	overall dimensions	
	mounting bracket)		
Platform Dimensions	60 in W x 32.5 in D x 2 in H / 1.524 mm W x 826 mm D x 51 mm H		
Ramp Dimensions	7.5 in W x 32.5 in D x 2 in H / 191 mm W x 826 mm D x 51 mm H		
Optional Guide	24 in W x .1 in D x 3.3 in H / 610 mm W x 3 mm D x 83 mm H		
Rail Dimensions	(2 short rails per scale) and 27 in W x .1 in D x 3.3 in H /		
	686 mm W x 3 mm D x 83 mm H (2 long rails per scale)		
Net Weight	115 lb / 52 kg (including	130 lb /	
	indicator and mounting bracket)	59 kg	
Shipping Weight	200 lb / 91 kg (including pallet)	215 lb / 98 kg (including pallet)	
Shipping	65 in L x 41 in W x 14 in H / 165 cm L x 104 cm W x 36 cm H		
Dimensions	USA		
Country of Origin			
UPC Code	809161305005	809161305104	

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



<sup>\*</sup> 203 E. Daugherty, Webb City, MO 64870 USA Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

Mobile: m.detecto.com

### **SOLD BY:**