

# VETERINARY PRODUCTS

Small animal and livestock scales, medical carts, waste receptacles, and glove box holders for veterinary practices



VET550



8440



VET75



**solace**  
IN-FLOOR ANALYSIS SCALES  
ID-4848S-855RMP



\* MOBILE AND PORTABLE \*  
LIVESTOCK SCALES

LSM1513



**RESCUE**   
medical carts  
RC33669BLU-L



GH2SS

PR Series



**DETECTO**<sup>®</sup>

# PRODUCTS FOR ALL SIZE ANIMALS, BIG OR SMALL

DETECTO manufactures a wide variety of scales, carts, glove box holders, and waste receptacles to outfit any veterinary practice or animal health facility. The most popular scale is the model VET550 which offers a 22 x 36 inch stainless steel portable platform, 550 lb / 250 kg capacity, and remote display for an economical price. DETECTO scales range in animal size from the small model 8440 with 44 lb x 1/2 oz capacity up to the big LSM1513 mobile livestock scale with 15,000 lb x 5 lb capacity.

## VET550

Medium-Size Animal Stainless Steel Scale with Remote Display



This economical, portable veterinary scale has a 22 in x 36 in stainless steel platform that is only two inches high. It's ideal for weighing dogs and cats and even goats and other smaller farm animals. The beautiful lightweight scale has built-in wheels and a guide handle making it ultra portable. Situate the remote display wherever you would like, near the platform or up on a countertop or wall mounted. It is powered by 4 AA alkaline batteries (not included) or the included AC adapter. 550 lb x 0.2 lb / 250 kg x 0.1 kg capacity.

## VET75

Precision, Portable Performance for Weighing Small Animals



DETECTO's compact VET75 digital veterinary scale is the right choice for weighing cats and other smaller size animals with its 75 lb x 0.02 lb, 704 oz x 0.5 oz, and 35 kg x 0.01 kg capacity. The all-in-one scale base, platform, and indicator is 20 in L x 14 in W x 2 in H / 50 cm L x 35 cm W x 5 cm H. The programmable Auto-Shutoff turns off scale when not in use to conserve batteries (adjustable up to 5 minutes). Veterinarians and vet techs will appreciate the TARE function, which can remove the weight of a container used to contain small animals while weighing.

## 8440

Small Animal Scale with Contoured Weighing Tray



Providing accurate and safe weighing of small animals such as cats, rabbits, gerbils, and hamsters, DETECTO's 8440 scale is ideal for veterinary clinics. The 8440's weighing tray may be removed to provide a flat weighing platform when needed and for easier transport. The 1-in / 25-mm high LCD offers pounds and kilograms readouts plus features a built-in clock. The scale's Weight Difference key provides the weight difference between the current weight on the scale and a saved weight in memory. Capacity is 44 lb x 1/2 oz / 20 kg x 10 g.

## RC33669BLU-L

RESCUE Steel Medical Carts



DETECTO's Rescue medical carts are sturdy and reliable... built to make sure you have the supplies you need where and when you need them for easy accessibility. They come in several colors: blue, red, white, and yellow. Carts have your choice of five or six drawers that can be secured with a key or EMG breakaway lock. The ABS countertop is easy to clean, and 5-inch wheels feature two parking locks and one steering lock. Optional accessories include IV poles, drawer dividers, glove dispensers, sharps container holders and waste bins.

## P-100R

Step-On Cans, the Rugged Waste Receptacle



Named for their rugged toughness, Step-On Cans are perfect for veterinary practices. Equipped with a heavy-duty foot pedal and linkage and a bag securing mechanism, they come in a variety of five different sizes and capacities (4, 6, 8, 12 and 25 gallon) and two different materials: stainless steel or a baked epoxy finish (available in red or white). The cans are flame retardant, have a self-closing lid for sanitation, and feature unitized construction.

## GH2SS


Wall-Mount Glove Box Holders



Conserve counter space with DETECTO's stainless steel or steel powder-coat painted universal glove box holders. The open design allows you to see what size or type of glove is being dispensed. The holders accommodate most glove boxes and are universal enough to be used for tissue holders where counter or shelf space is limited. GH models are available for holding one, two, or three glove boxes. Lifetime warranty.

## ID-4848S-855RMP

Solace In-Floor Flush-Mounted Scale

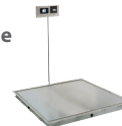

Solace series in-floor scales are flush-mounted platform scales that feature stainless steel trim decks and include a pit frame, trim ring, and interface cable to the included multi-color, touchscreen LCD indicator and printer. The flush surface eliminates tripping hazards and hallway accessibility issues and provides you with a scale that can weigh virtually any indoor-size animal. The Solace includes a recessed-mount indicator and ticket printer combination, high 1,000 lb x 0.2 lb / 450 kg x 0.1 kg capacity, and they come in a variety of platform sizes from 3 x 3 to 6 x 4 feet.

## LSM1513-MWC225P

Weight Wrangler Mobile Group Livestock Scale



Weight Wrangler livestock scales provide a mobile cattle scale that can weigh one cow or up to 15 all at once. The LSM series can be pulled behind a 3/4-ton or larger truck and provides a self-contained NTEP legal-for-trade weighing system that can be used on any firm surface with up to a 5% grade. Models come in 18-ft-long tandem axle or 13-ft-long single axle sizes. The scale includes a weight indicator, ticket printer, USB port, battery, self-adjusting electric brakes, and adjustable hitch height. Solar panels are optional.

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

**SOLD BY:**

# DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA  
Toll-Free: (800) 641-2008 E-mail: [detecto@cardet.com](mailto:detecto@cardet.com)  
[www.Detecto.com](http://www.Detecto.com)