

DETECTO®

WPS12UT

DIGITAL SCALE WITH UTILITY BOWL



WPS12UT

WEIGH




Chux Pads



Wet Diapers



Organs

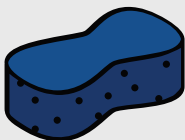

Lap Sponges

Submersible Underwater for Hygienic Cleaning

- Accurately measure fluids or blood loss during clinical procedures
- High-quality stainless steel construction
- Rechargeable battery pack included
- Accurate 1 gram division size
- Six different capacities, including 192 oz x 0.05 oz and 5,500 g x 1 g

COMPLETELY SUBMERSIBLE UNDER WATER

FOR EASY
CLEANING



IDEAL
FOR
QBL
WEIGHT

Excellent for weighing diaper fluids or blood loss absorbed by chuxs used in the OR when mothers hemorrhage during labor.

Completely-sealed design allows the scale and bowl to be submersible underwater for sanitary washdown environments

DETECTO's model WPS12UT IP67-rated washdown digital scale with utility bowl can be submersed underwater for easy clean up that is necessary in clinical industry applications. This scale is ideal for weighing used chux underpads, wet diapers (for both adults and children), lap sponges, and organs.

- 4.6 in / 11.7 cm deep bowl aids sample retention and allows easy clean-up
- Tare button removes the dry diaper or chux weight and measures fluid output of the soiled diaper or bloody chux
- Built-in rechargeable battery pack and AC adapter included
- 12-lb capacity with 0.1-oz division size



- Attractive backlit blue LCD readouts with convenient display viewing angle
- Stainless steel scale and bowl for sterile environments



The stainless steel scale and bowl construction prevents corrosion and keeps water from reaching the more delicate working parts of the scale electronics.



WPS12UT

DIGITAL SCALE WITH UTILITY BOWL



**WEIGH, MEASURE, AND
QUANTIFY BLOOD LOSS
FOR OPTIMAL QBL**



Whether the fluid loss is collected by lap sponges, diapers, or underpads, DETECTO's WPS12UT digital scale will weigh it accurately and precisely in a second.

- Raised Buttons Allow for Tactile Functionality and Easy Operation
- Convenient Upward-Facing Weight Display Viewing Angle
- Tare Function to Remove the Weight of a Dry Chux or Diaper
- Instantly Displays Weight in Both Decimals and Fractions
- 1.2-in/30-mm High Backlit LCD is Easily Readable in Dimly-Lit Areas
- Toggle Through Six Different Weight Capacity Display Options
- Convenient 80-Hour Rechargeable Battery Pack Allows the Scale to be Moved Wherever it's Needed Most



RECHARGEABLE BATTERY & WATERTIGHT PLUG

The NSF-certified stainless steel scale offers a dual power source of rechargeable battery power or 100-240V AC adapter (included) with a watertight AC plug on the side. The automatic shutoff feature helps conserve battery life.

DETECTO

WEIGH




Chux Pads



Wet Diapers



Organs

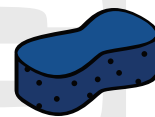

Lap Sponges

LARGE 1-INCH-HIGH LCD DIGITS

GRAMS FRACTIONAL OUNCES AND MORE

The WPS12UT includes the submersible WPS12 Mariner® scale, roomy stainless steel bowl, AC adapter, quick set-up instruction sheet, and full operation manual.



NSF Certified

The NSF-certified WPS12UT saves time and helps reduce the risk of bacteria harboring in the nooks and crannies of the scale by providing complete submersion underwater when washing for sanitary environments.



	WPS12UT
Capacity/Resolution	12 lb x 0.1 oz / 12 lb x 0.005 lb / 12 lb x 1/8 oz / 192 oz x 0.05 oz / 192 oz x 1/8 oz / 5,500 g x 1 g
Overall Dimensions (Scale and Bowl)	13.6 in W x 13.6 in D x 6.9 in H / 34.5 cm W x 34.5 cm D x 17.5 cm H
Bowl Dimensions	13.6 in W x 13.6 in D x 4.6 in H / 34.5 cm W x 34.5 cm D x 11.7 cm H
Construction	Stainless steel (bowl and scale enclosure)
Display	1.2-in / 30-mm high backlit blue LCD, 5 digits
Digits	1-in / 25-mm high black LCD
Keypad Buttons	ON/OFF, UNITS/Backlight, SET, TARE
Power	AC Adapter (with US plug included): 100-240 VAC input, 9V @ 1A output Battery (included): Ni-MH rechargeable 4.8 VDC 1800 mA (Note: a plug adapter is required for international use)
Battery Life	80 hours of continuous use with a fully recharged battery pack
Feet	6 rubberized feet for traction
Tare	100% of full scale capacity
Zero	Established on power up routine and maintained by auto-zero circuitry
Enclosure Rating	IP67 submersible
Certification	NSF Approved
Operating Temperature	40° to 105°F (5° to 40°C)
Storage Temperature	40° to 140°F (5° to 60°C). Please Note: The WPS12UT is not commercial dishwasher safe.
Net Weight	4.6 lb / 2.1 kg
Shipping Weight	7.6 lb / 3.4 kg
UPC Code	809161317404

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

SOLD BY:



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com