

DETECTO®

apex®

CLINICAL SCALES

EYE-LEVEL DIGITAL CLINICAL SCALES

- 600 lb x 0.1 lb / 270 kg x 0.05 kg capacity
- Sonar touchless or inline mechanical height rod
- Extra-wide and flat platform (17 in W x 17 in D x 2.75 in H)
- Wi-Fi/Bluetooth models available
- Units locking into kilograms or pounds



Model APEX-SH

Model APEX



**WEIGHT + HEIGHT
INSTANT READOUTS
WITHOUT PUSHING
A SINGLE BUTTON**

Sonar touchless height rods available on some apex® models provide immediate digital readouts for optimal speed, accuracy, and hygiene

THE APEX[®]

OF PATIENT
MEASUREMENT
EFFICIENCY
+ STYLISH
DESIGN
IS NOW AVAILABLE!

STATE-OF-THE-ART SONAR HEIGHT RODS

DETECTO's unique touchless sonar height rod instantly locks on the patient's height and measures nearly any patient size. The apex's sonar height rod measures from 24 inches (2 ft) up to 86 inches (7'2") by 0.2 inch increments and metric readouts of 60 cm up to 218 cm by 0.5 cm increments (model APEX-SH). Simply step on the scale and your height will be captured and displayed on the LCD within a second.




DETECTO's apex[®]
scales measure weight
and height efficiently
on patients of all
ages and sizes.

The apex's inline mechanical height rod models measure from 24 inches (2 ft) up to 86 inches (7'2") by 0.125 inch increments and metric readings of 60 cm up to 218 cm by 0.1 cm increments (model APEX). To calculate BMI, simply use the scroll arrows to enter the height measurement and then press the BMI key.

- Eye-level weight display convenience
- Sturdy column construction
- AC adapter included
- 10 different standard apex[®] models available for you to choose from
- OneWeigh units locking into lbs or kgs



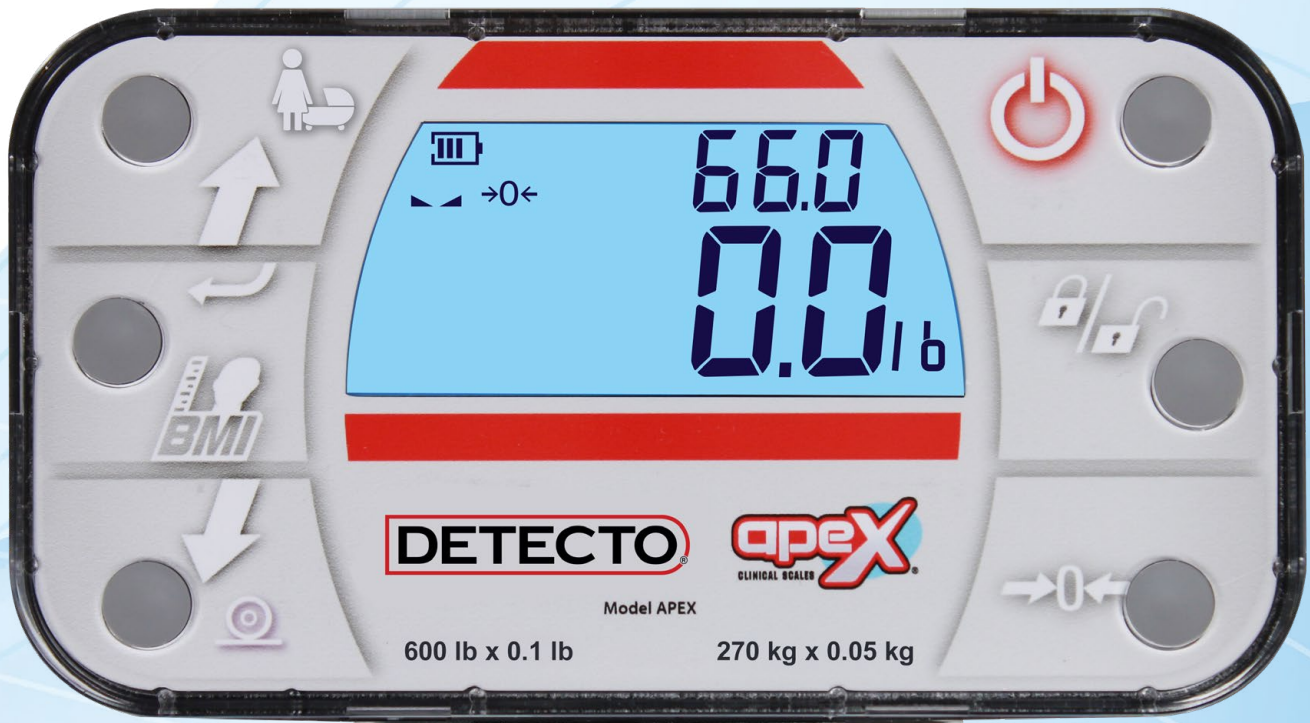

HIGH-TECH WEIGHT INDICATOR

The apex® features a high-tech, smart-phone-style indicator for easy utilization with 6 simple buttons and a bold LCD display for large readouts.

- 0.75-inch-high, clinical-blue LCD weight readouts
- Body Mass Index calculation
- Weight and height displayed onscreen simultaneously
- Automatically locks on weight temporarily for a few seconds in order to easily view and record measurement
- Battery power level indication
- Mother/Baby function factors the weight of infants and toddlers held by an adult



Wi-Fi/Bluetooth models available

Dual integral wheels provide for easy transport on all apex® models.



VERSATILE CONNECTIVITY:

- 1 micro USB-B port
- 1 RS232 serial port
- HL7 IEEE 11073 compliant (standard)
- Wi-Fi/Bluetooth models available for wireless EMR/EHR
- Welch Allyn® device connectivity available: CVSM, CSM, and LXI (see models chart on back page)




POWER



apex® models come standard with the AC adapter included (see models chart on back).

An integral bubble level helps locate the most level surface area on your floor.

The lift-off apex® platform cover may be removed for wipe-down cleaning or to replace the 12 AA batteries. The scale's structural base integrity is manufactured rock-solid for optimal patient safety and stability.

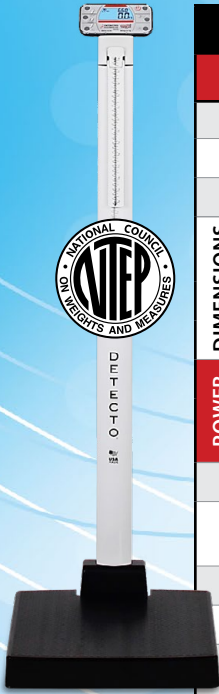

apex® EYE-LEVEL DIGITAL SCALES WITH MECHANICAL HEIGHT RODS					
	Model	APEX-AC	APEX-C-AC	APEX-UWA-AC	APEX-LXI-AC
	Capacity	600 lb x 0.1 lb / 270 kg x 0.05 kg			
	Height Rod	Mechanical inline height rod			
	Height Rod Measures	24" (2') - 86" (7'2") x 0.125 in / 60 cm - 218 cm x 0.1 cm			
DIMENSIONS	Overall Dimensions	17 in W x 19.5 in D x 53.3 in H / 43 cm W x 49.5 cm D x 135.3 cm H			
	Platform Size	17 in W x 17 in D x 2.75 in H / 43.2 cm W x 43.2 cm D x 7 cm H			
	LCD	.75 in / 1.9 cm H LCD digits			
POWER	Batteries	12 AA batteries (not included)			
	AC Adapter (included)	•	•	•	•
	Ports	1 RS232 serial port and 1 micro USB-B port			
	Connectivity		Wi-Fi / Bluetooth	Welch Allyn CVSM/CSM	Welch Allyn LXI
	Net Weight	31 lb / 14 kg			
	Shipping Weight	40 lb / 18 kg			
	Shipping Dimensions	58 in x 20 in x 6 in / 147 cm x 51 cm x 15 cm			
	UPC Code	809161197907	809161322606	809161198102	809161198201

apex® EYE-LEVEL DIGITAL SCALES WITH SONAR HEIGHT RODS					
	Model	APEX-SH-AC	APEX-SH-C-AC	APEX-SH-UWA-AC	APEX-SH-LXI-AC
	Capacity	600 lb x 0.1 lb / 270 kg x 0.05 kg			
	Height Rod	Touchless sonar height rod			
	Height Rod Measures	24" (2') - 86" (7'2") x 0.2 in / 60 cm - 218 cm x 0.5 cm			
DIMENSIONS	Overall Dimensions	17 in W x 19.5 in D x 89.5 in H / 43 cm W x 49.5 cm D x 227.3 cm H			
	Platform Size	17 in W x 17 in D x 2.75 in H / 43.2 cm W x 43.2 cm D x 7 cm H			
	LCD	.75 in / 1.9 cm H LCD digits			
POWER	Batteries	12 AA batteries (not included)			
	AC Adapter (included)	•	•	•	•
	Ports	1 RS232 serial port and 1 micro USB-B port			
	Connectivity		Wi-Fi / Bluetooth	Welch Allyn CVSM/CSM	Welch Allyn LXI
	Net Weight	32 lb / 14 kg			
	Shipping Weight	41 lb / 19 kg			
	Shipping Dimensions	58 in x 20 in x 6 in / 147 cm x 51 cm x 15 cm			
	UPC Code	809161197907	809161322804	809161198102	809161198201



OIML Class III Compliant: These models are also available as OIML Class III compliant by adding -KGEU suffix to the model number when ordering. Capacity: 290 kg x 0.1 kg.

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.

DETECTO®

102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com

SOLD BY: