

DETECTO®



SOLO®
DIGITAL CLINICAL SCALES

DETECTO'S SOLO® SCALE STANDS ALONE IN ECONOMIC CLINICAL WEIGHING

- 550 lb x 0.2 lb / 250 kg x 0.1 kg Capacity
- Weight, Height, and Body Mass Index
- 14 in W x 15 in D x 2.5 in H Platform Size
- 6 AA Batteries (included) or
Optional AC Adapter
- Inline Mechanical Height Rod
- USB Port for EMR/EHR Connectivity
- Mother/Baby Mode Factors the
Weight of Infants Held by an Adult



**UNMATCHED
ECONOMICAL
WEIGHING
QUALITY**



*The Name Brand Preferred
by Clinical Professionals
Around the World*



Eye-level height display for ease of use (52 in / 132 cm H)

Inline height rod measures:
30 in to 84.5 in /
76 cm to 214.5 cm

1-inch/25-mm high LCD,
clearly-visible display



OneWeigh units locking into Pounds or Kilograms avoids dosage errors based on incorrect units.

Rock-solid, stable column – aluminum white powder coated

Two wheels for quick and easy portable transport

Rubber mat for patient assurance and stability

DETECTO



Stylish Indicator

The solo® indicator features 8 simple buttons that are easily identifiable in any language.



Battery Power

No power outlet available, no problem! The solo® comes with 6 AA batteries for weighing anywhere.




EMR/EHR-Ready

The scale's USB port allows you to output measurement data to EMR/EHR software.




Body Mass Index

Simply step on the scale to weigh and the scroll arrows let you rapidly enter your height to factor BMI.

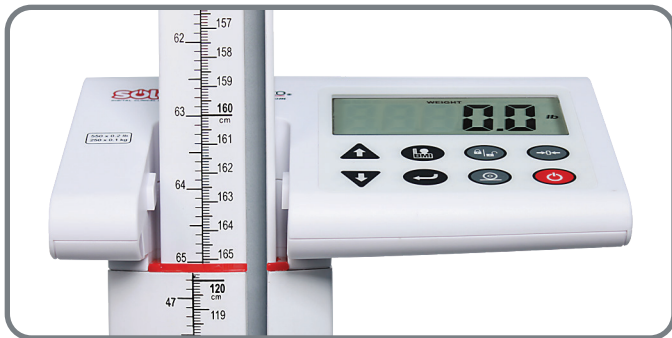

Inline Height Rod

The mechanical height rod is situated to the side of the LCD, so it never interferes with viewing measurements.

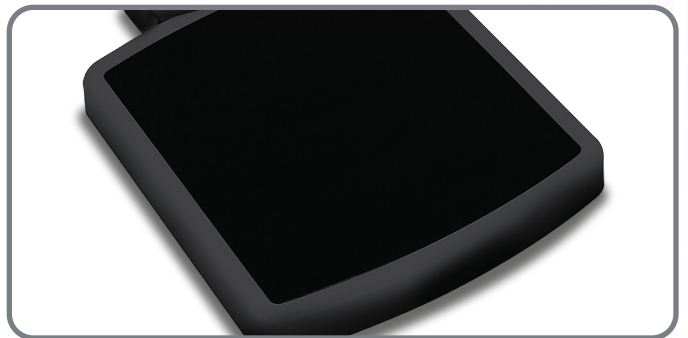

Wide, Flat Platform

The oversized 14 x 15 inch platform with rubber mat is only 2.4 inches high making it easy for patients to step on/off.



AC Adapter

The AC port is located in the base, so no pesky wires dangle from up on the indicator (AC Adapter optional: model PD-AC).



Portable, Wheels

Two integral wheels allow for ease of movement whenever you need to transport the solo®.



Model	solo
Capacity	550 lb x 0.2 lb / 250 kg x 0.1 kg
Platform Size	14 in W x 15 in D x 2.4 in H / 36 cm W x 38 cm D x 6 cm H
Weight Units	Pounds or Kilograms (selectable upon initialization at start-up; OneWeigh units locking to avoid dosage errors)
Power	Six (6) AA size Alkaline batteries (included) OR an optional medical device 9V AC/DC wall plug-in adapter (Detecto part number PD-AC). The AC adapter comes with model SOLO-AC.
Display	Five-digit, seven-segment, 1.0-inch (25-mm) high LCD digits
Operation Temperature	50 to 104 °F / +10 to +40 °C
Function Keys	ON/OFF, ZERO, BMI HEIGHT, ENTER, LOCK/RELEASE, PRINT, UP ARROW/Mother and Baby, DOWN ARROW Interface: USB (ASCII Protocol)
Annunciators	ZERO, LOCK, STA, HOLD, LB, KG, CM, FT/IN, HEIGHT, WEIGHT, BMI, AM, PM
Scale Dimensions	14.3 in W x 19.5 in D x 52.0 in H / 36.3 cm W x 49.5 cm D x 132.1 cm H
Interface	USB (ASCII Protocol)
Height Rod	Mechanical included standard, inline
Height Rod Units	Inches or Centimeters
Height Rod Range	30 in to 84.5 in / 76 cm to 214.5 cm
Manual Input Height	2 ft 0.0 in – 7 ft 10.5 in by 0.5 in / 60.0 cm - 240.0 cm by 1cm
Minimum Weighing	11 lb / 5 kg
Wheels	(2) wheels for ease of movement
Selectable Auto Shut-Off	Up to 9 Minutes
Construction	Column and height rod – aluminum alloy white powder coated; indicator enclosure, height rod arm, platform and base – white ABS; mat – black PVC
Product Weight	17.7 lb / 7.7 kg
Shipping Weight	25 lb / 11.3 kg
Package Dimensions	30 in W x 20 in D x 11 in H / 76 cm W x 51 cm D x 28 cm H
UPC Code	809161197709



DETECTO's SOLO-RI is also available with an eye-level column as model SOLO

Model SOLO

Model SOLO-RI



102 East Daugherty St., Webb City, MO 64870 USA
Toll-Free: (800) 641-2008 E-mail: detecto@cardet.com
www.Detecto.com

SOLD BY: